



THE PT ELEVATOR™

PT Elevator Owner's Safety and User Guide

SAFETY

FOR CHILDREN 12 AND UNDER, USE UNDER SUPERVISION ONLY.

The intended use of this product is to raise and lower a person to and from the floor.

NOTE: This owner's manual is compiled from the latest specifications and product information available at the time of publication. We reserve the right to make changes as they become necessary. Any changes to our products may cause slight variations between the illustrations and explanations in this manual and the product you have purchased. The latest/current version of this manual is available on our website.

Read and follow all instructions, warnings, and notes in this manual before attempting to operate your chair for the first time. Your safety depends upon you, as well as your provider, caretaker, or healthcare professional in using good judgment.


KEEP HANDS AND FEET AWAY FROM MOVING PARTS AND EDGE OF SEAT (PINCH POINTS).



Do not, under any circumstances, tamper with the batteries. You should use the charger supplied with the chair. It is the safest, most efficient tool to charge the batteries. We do not recommend using other types of chargers.

Charging the Chair

- > Connect the battery charger round output plug to the charging socket on the side panel of chair
- > Connect the charger to a wall outlet by means of the main plug.

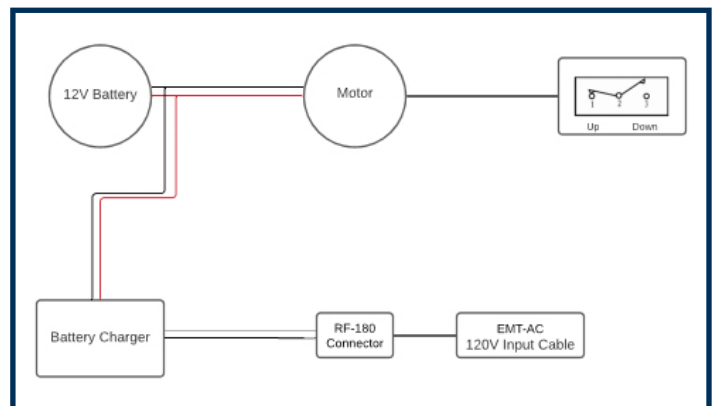



USAGE INSTRUCTIONS

FOR CHILDREN 12 AND UNDER, USE UNDER SUPERVISION ONLY.

The intended use of this product is to raise and lower a person to and from the floor.

Read and follow all instructions, warnings, and notes in this manual before attempting to operate your chair for the first time. Your safety depends upon you, as well as your provider, caretaker, or healthcare professional in using good judgment.

Lowering to the Floor

1. Press and hold toggle switch to raise or lower seat to required height to sit easily on the seat. NOTE: Be sure to sit completely on seat, with back against the backrest.
2. Once comfortably seated, press and hold toggle switch to lower chair to floor level.

Rising from the Floor

1. When seated completely and comfortably on the chair, press and hold the toggle switch to raise the chair to required height to easily stand up from chair.

Moving the Chair

To engage wheels:

- > Raise or semi-raise chair
- > Tilt chair backwards
- > Roll to desired area
- > Lower chair back to a level position, fully on base



Designed, patented, and made in the USA.

For more information, please visit www.rehabchairlift.com or email info@rehabchairlift.com.